

SCIENTIFIC REPORTS

OPEN

Mycobacterial genomic DNA from used Xpert MTB/RIF cartridges can be utilised for accurate second-line genotypic drug susceptibility testing and spoligotyping

Rouxjeane Venter¹, Brigitta Derendinger¹, Margaretha de Vos¹, Samantha Pillay^{1,2}, Tanya Dolby², John Simpson², Natasha Kitchin¹, Ashley Ruiters², Paul D. van Helden¹, Robin M. Warren¹ & Grant Theron¹

Xpert MTB/RIF (Xpert) is a widely-used test for tuberculosis (TB) and rifampicin-resistance. Second-line drug susceptibility testing (DST), which is recommended by policymakers, typically requires additional specimen collection that delays effective treatment initiation. We examined whether cartridge extract (CE) from used Xpert TB-positive cartridges was, without downstream DNA extraction or purification, suitable for both genotypic DST (MTBDR_{plus}, MTBDR_{sl}), which may permit patients to rapidly receive a XDR-TB diagnosis from a single specimen, and spoligotyping, which could facilitate routine genotyping. To determine the limit-of-detection and diagnostic accuracy, CEs from dilution series of drug-susceptible and -resistant bacilli were tested (MTBDR_{plus}, MTBDR_{sl}). Xpert TB-positive patient sputa CEs (n = 85) were tested (56 Xpert-rifampicin-susceptible, MTBDR_{plus} and MTBDR_{sl}; 29 Xpert-rifampicin-resistant, MTBDR_{sl}). Spoligotyping was done on CEs from dilution series and patient sputa (n = 10). MTBDR_{plus} had high non-valid result rates. MTBDR_{sl} on CEs from dilutions $\geq 10^3$ CFU/ml ($C_T \leq 24$, > "low" Xpert semiquantitation category) was accurate, had low indeterminate rates and, on CE from sputa, highly concordant with MTBDR_{sl} isolate results. CE spoligotyping results from dilutions $\geq 10^3$ CFU/ml and sputa were correct. MTBDR_{sl} and spoligotyping on CE are thus highly feasible. These findings reduce the need for additional specimen collection and culture, for which capacity is limited in high-burden countries, and have implications for diagnostic laboratories and TB molecular epidemiology.

Of the 10.4 million individuals with active tuberculosis (TB) in 2015, 580 000 were rifampicin (RIF) resistant or multidrug-resistant (MDR), defined as resistance to isoniazid (INH) and RIF¹. Only ~20% of MDR-TB cases were diagnosed and started on treatment, and only half started on treatment were cured¹. Extensively drug-resistant (XDR)-TB, which is MDR with resistance to a fluoroquinolone (FQ) and a second-line injectable drug (SLID) comprises 10% of MDR-TB cases, yet is even more underdiagnosed than MDR-TB, very costly to treat, and represents an emerging public health emergency²⁻⁶.

Xpert MTB/RIF (Xpert) (Cepheid, United States) is a Food and Drug Administration and World Health Organization (WHO)-endorsed nucleic acid amplification test (NAAT) that rapidly detects *Mycobacterium tuberculosis* complex-DNA and RIF-resistance directly from sputa⁷⁻⁹. Over 25 million Xpert MTB/RIF cartridges have been consumed and over 30 000 test modules are installed worldwide¹⁰. The WHO and several national programmes recommend that if Xpert detects resistance, an additional sputum is collected for further drug

¹DST/NRF Centre of Excellence for Biomedical Tuberculosis Research, SA MRC Centre for Tuberculosis Research, Division of Molecular Biology and Human Genetics, Faculty of Medicine and Health Sciences, Stellenbosch University, Cape Town, South Africa. ²National Health Laboratory Services, Cape Town, South Africa. Rouxjeane Venter and Brigitta Derendinger contributed equally to this work. Correspondence and requests for materials should be addressed to G.T. (email: gtheron@sun.ac.za)

susceptibility testing (DST) using line probe assays (LPAs), such as MTBDR*plus* (RIF and INH) and MTBDR*sl* (FQs and SLIDs), or phenotypic testing^{1,9,11,12}.

Patients, however, often do not rapidly return to the clinic to give another sputum or receive DST results. For example, a study in South Africa found that, even after MTBDR*plus* roll-out, time-to-treatment since initial diagnosis was unacceptably long (~55 days), and that this was partly due to challenges with patient loss-to-follow-up¹³. Furthermore, many patients do not produce sufficient sputum of adequate quality, especially in settings with high rates of HIV^{14–17}.

MTBDR*plus* and MTBDR*sl* have suboptimal sensitivity on specimens, and culture is often required prior to DNA extraction and further genotypic testing. Not only can this cause diagnostic delay, but many high burden countries lack the necessary biosafety and laboratory infrastructure for mycobacterial culture and DNA extraction^{18–21}. Furthermore, culture can result in the loss of potentially clinically-meaningful resistance²². There is hence a need to reduce delays in the diagnosis of drug-resistant TB and use rapid methods that minimise reliance on culture through the direct testing of specimens²³.

Poor adherence to diagnostic algorithms using MTBDR*plus* and MTBDR*sl* has been reported^{5,24,25}. For example, in South Africa, 34% of Xpert RIF-resistant patients failed to receive MTBDR*plus* and, of those confirmed to have MDR-TB, 28% did not receive second-line DST with MTBDR*sl* – despite both LPAs being mandated by the national programme²¹. Novel approaches to reduce this gap in the TB care cascade, which is worsened by the requirement for extra patient visits and additional specimen collection, is a major research priority^{26,27}. If TB-testing and first- and second-line DST were possible on the first available specimen, fewer patients would potentially be lost and patients could be diagnosed earlier. This could result in earlier effective treatment initiation, fewer patient- and health systems-costs, and better long-term clinical outcomes.

We therefore conducted a proof-of-concept evaluation on whether *M. tuberculosis*-complex genomic DNA in the PCR-reaction mix from used Xpert cartridges (cartridge extract; CE) – that would otherwise be discarded – was detectable in an accurate manner using MTBDR*plus* and MTBDR*sl*. The feasibility of genotyping on CE by spoligotyping was also tested as this would potentially be useful for research laboratories and programmes seeking to implement routine strain surveillance. We explored the feasibility of Sanger sequencing on CE, as this may be useful for additional genotypic DST. Critically, we evaluated CE for all tests without additional downstream DNA extraction or purification, as not only would extraction require equipment not readily available in routine diagnostic laboratories in high burden settings, but it would complicate laboratory workflows and reduce the attractiveness of our approach. If the CE approach was feasible, it would mean that many laboratories would already have instrumentation available for mycobacterial genomic DNA extraction in the form of GeneXpert¹⁰ and not need to procure new equipment.

Material and Methods

Ethics statement. Methods and protocols were carried out in accordance with relevant guidelines and regulations. The study was approved by the Health Research Ethics Committee of Stellenbosch University (N09/11/296) and the City of Cape Town (#10570). Permission was granted to use anonymised residual specimens collected as part of routine diagnostic practice and thus informed consent was waived.

Xpert MTB/RIF on dilution series of drug-susceptible- and drug-resistant bacilli. A triplicate ten-fold dilution series was made using phenotypically-confirmed drug-susceptible (DS)-TB, MDR-TB and XDR-TB clinical isolates (0–10⁶ CFU/ml) in phosphate buffer (33 mM Na₂HPO₄, 33 mM KH₂PO₄; pH 6.8) with 0.025% Tween80 (Sigma-Aldrich, United States). Colony counts were done by plating on 7H11 Middlebrook agar (BD Biosciences, United States). Dilutions containing bacilli (1 ml aliquots) were tested by Xpert (54 in total: six dilutions ranging from 10¹–10⁶ CFU/ml in triplicate for three strains and hence 18 dilutions each for the DS, MDR, and XDR strains) as well as 0 CFU/ml controls in triplicate, according to the manufacturer's instructions⁹. Used cartridges were stored at 4 °C prior to CE extraction within 24 h and freezing of the CE at –20 °C.

Xpert MTB/RIF on clinical specimens. Used Xpert-TB-positive cartridges done on sputa from people with symptoms suggestive of TB tested as part of the South African national TB diagnostic algorithm were collected between February 2016 and November 2016 from the National Health Laboratory Services (NHLS), a South African National Accreditation System-accredited, quality-assured diagnostics laboratory in Cape Town, South Africa¹¹. Cartridges were stored at 4 °C prior to CE extraction within 5 days. Fifty-six Xpert TB-positive, RIF-susceptible cartridges and 29 Xpert-TB-positive RIF-resistant cartridges were collected. When the NHLS did a MGIT 960 liquid culture on sputum from RIF-resistant patients, we collected the isolate [20/29 (69%) had available isolates]. Isolates were not available from Xpert TB-positive, RIF-susceptible specimens as culture is not routinely done in these patients^{11,28}.

Recovery of mycobacterial genomic DNA from used Xpert MTB/RIF cartridges. The transparent diamond-shaped reaction chamber on the back of the cartridge was punctured with a sterile fixed-needle insulin syringe (1 ml; 29 G) (Fig. 1) in a biosafety level 2 cabinet. The full CE volume, typically ~15 µl, was withdrawn and stored in sterile, safe-lock micro-centrifuge tubes at –20 °C prior to analysis. Each cartridge and the surrounding surface was wiped down thoroughly with 1% sodium hypochlorite and 70% EtOH before and after extraction and UV sterilization was done after each batch of extraction. Used needles were discarded in a sharps container containing 1% sodium hypochlorite. Before and after each cartridge extraction session, hood surface area was decontaminated with sodium hypochlorite and EtOH and UV sterilised. No DNA extraction or purification steps were done on CE.

Line probe assays on cartridge extract. MTBDR*plus* and MTBDR*sl* (both version 2.0) were done according to the manufacturer's instructions^{29,30} except for Xpert TB-positive, RIF-susceptible clinical specimens CE (n = 56), 7.5 µl CE was used as input volume into MTBDR*plus* and MTBDR*sl*. For the Xpert TB-positive,



Figure 1. Cartridge extract extraction procedure. (a) The arrow indicates the diamond-shaped reaction chamber where the PCR amplification takes place and contains cartridge extract with mycobacterial genomic DNA. The needle is placed at the top of the diamond and the film is slowly and carefully pierced. (b) The needle is then slowly inserted deeper into the pocket and cartridge extract mix drawn out without piercing the other side.

RIF-resistant clinical specimen CEs (n = 29) and the dilution series, 5 µl (the recommended input volume) CE was used in order to have enough CE remaining for Sanger sequencing. MTBDR*plus* and MTBDR*sl* results were reported as susceptible or resistant (RIF and INH for MTBDR*plus*; FQ and SLID for MTBDR*sl*), indeterminate [*M. tuberculosis* complex DNA-positive (reported by the test as TUB-positive) but no gene loci control bands] or TUB-band negative. LPA strips were read by two independent, experienced readers blinded to each other's calls and Xpert results (and, for dilution series, the strain used).

Spoligotyping on cartridge extract. Spoligotyping was done as described^{31,32} on 2 µl CE from the MDR-TB dilution series. A set of Xpert TB-positive, RIF-susceptible cartridges (n = 10) done on specimens and separate from those used for genotypic DST on CE were collected with paired culture isolates from an ongoing research study. To determine whether the correct spoligotype was obtained from CE, crude DNA extracted through heat inactivation from the corresponding culture isolates was spoligotyped. SITVIT was used to identify strain families³³.

Targeted Sanger sequencing on cartridge extract. For dilution series, PCR clean-up and Sanger sequencing on 5 µl CE was done by the Stellenbosch University Central Analytical Facility using primers overlapping with LPA-binding sites (Supplementary Table 1). The *gyrA* and *rrs* regions in the DS-TB and XDR-TB strains were sequenced.

Data availability. The datasets generated during and/or analysed during the current study are available from the corresponding author on reasonable request.

Results

Patient characteristics. A summary of the patient demographic and clinical data is in Table 1. For Xpert TB-positive, RIF-susceptible patients the median age (IQR) was 40 (31–49) years and for RIF-resistant specimens was 35 (23–42) years. 37/55 (67%) of RIF-susceptible patients and 12/20 (60%) of RIF-resistant patients were male.

Feasibility and diagnostic accuracy of MTBDR*plus* and MTBDR*sl* on dilution series Xpert TB-positive cartridge extract. Xpert detected *M. tuberculosis*-complex DNA in all dilutions $\geq 10^2$ CFU/ml and correctly identified RIF-susceptibility and -resistance (Fig. 2). MTBDR*plus* showed poor overall sensitivity for *M. tuberculosis*-complex DNA [22% (12/54) TUB-band-positive] in CE from Xpert TB-positive cartridges. MTBDR*plus* had high rates of non-actionable (TUB-band negative or TUB-band positive but indeterminate) and false RIF-heteroresistant results (Figs 2 and 3).

In contrast, MTBDR*sl* on CE had high sensitivity and specificity [87% (47/54) and 100% (9/9) respectively] for *M. tuberculosis*-complex DNA and a limit of detection of 10^3 CFU/ml. Susceptibility and resistance to FQs and SLIDs were correctly detected for all strains $\geq 10^3$ CFU/ml, corresponding to $C_T \leq 24$ (the higher C_T range of the Xpert “low” semiquantitation category) in all but one sample (one replicate of the MDR-TB strain was indeterminate for FQs; Fig. 3). Once non-actionable results were excluded, overall sensitivities and specificities of 87% (13/15) and 96% (25/26) for FQ-resistance and 94% (15/16) and 97% (30/31) for SLID-resistance, respectively were obtained. When the threshold of $\geq 10^3$ CFU/ml ($C_T \leq 24$) was applied, the sensitivity and specificity were both 100% (12/12 and 23/23, respectively) for FQs and for SLIDs (12/12 and 24/24, respectively).

| Patient Characteristics | Xpert TB-positive | |
|-----------------------------------|---------------------------------------|-------------------------------------|
| | Xpert rifampicin-susceptible (n = 56) | Xpert rifampicin-resistant (n = 29) |
| Age, median (IQR) | 40 (30–49) | 35 (23–42; p = 0.086) |
| Male gender (%) | 37/55 (67)* | 12/20 (60)* |
| Smear-positivity (%) | 37/50 (74)* | 6/16 (38)* |
| Culture-positivity (%) | Not done | 19/21 (90)* |
| TTP, median (IQR) | N/A | 10 (8–20) |
| Xpert C _T , median IQR | 17.9 (16.3–22.1) | 20.5 (16.9–24.8) |

Table 1. Patient demographic and clinical data. *Missing data: Gender (n = 1 for RIF-susceptible, n = 9 for RIF-resistant); Smear status (n = 6 for Xpert RIF susceptible, n = 13 for Xpert RIF-resistant); Culture results (n = 8 for RIF-resistant results). Abbreviations: Xpert - Xpert MTB/RIF; IQR - interquartile range; TTP - time-to-positivity; C_T - cycle threshold values.

Diagnostic accuracy of MTBDR_{plus} and MTBDR_{sl} on clinical specimen Xpert TB-positive cartridge extract. *Xpert MTB/RIF rifampicin-susceptible specimens.* As with the dilution series, MTBDR_{plus} had high rates of indeterminate and false-resistance results on clinical specimen CE (Table 2). However, most MTBDR_{sl} results from Xpert TB-positive, RIF-susceptible clinical specimen CE were valid (TUB-positive, not indeterminate, and no false-susceptible or -resistant results): 53/56 (95%) for FQ (two TUB-band negative, one indeterminate) and 51/56 (91%) for SLID (two TUB-band negative, three indeterminate). The few CEs that yielded indeterminate MTBDR_{sl} results corresponded to “low” or “very low” Xpert semiquantitation levels (C_T > 24). The median (IQR) C_T of indeterminate (26.3, 24.4–26.7) vs. determinate (17.62, 15.6–20.6) MTBDR_{sl} results differed (p < 0.001), indicating that indeterminate results are likely a function of low DNA concentrations in CE. There was not enough CE volume remaining or a matching clinical isolate for confirmatory testing from the three MTBDR_{sl}-detected SLIDs resistant patients.

Xpert MTB/RIF rifampicin-resistant specimens. MTBDR_{sl} on Xpert TB-positive, RIF-resistant CE had 24/29 (83%) valid results. For FQs, 14/24 (58%) were susceptible and 10/24 (42%) were resistant. For SLIDs, 15/24 (63%) were susceptible and 9/24 (37%) resistant. The five non-valid results were TUB-band-negative [2/29 (7%)] or indeterminate for both FQs and SLIDs [3/29 (10%); Table 2]. All CEs corresponding to the higher C_T ranges of the Xpert “low” semiquantitation category (C_T ≤ 24) had valid results, whereas those that had indeterminate or TUB band-negative results corresponded to the lower semiquantitation levels (C_T > 25.0). The median (IQR) C_T of indeterminate (29.1, 26.5–31.1) vs. determinate (20.5, 16.–23.2) results differed significantly (p < 0.001).

MTBDR_{plus} and MTBDR_{sl} performance on Xpert MTB/RIF cartridge extract by smear status. MTBDR_{plus} had high non-valid result rates irrespective of smear status. However, MTBDR_{sl} on CE from smear-negative sputums had significantly higher rates of non-actionable results [5/23 (22%) vs. 1/43 (2%) for FQ, p = 0.01; 6/23 (23%) vs. 2/43 (5%) for SLIDs, p = 0.01] compared to smear-positive patients (Supplementary Table 2).

Concordance of MTBDR_{sl} results on cartridge extract and culture isolates. Of the 29 Xpert TB-positive, RIF-resistant patients, 20 (69%) matched culture isolates were collected while the remaining nine had negative or contaminated cultures. The CEs and isolates showed 18/20 (90%) matching MTBDR_{sl} FQ results and 17/20 (84%) matching SLID results. There were 2/20 (10%) discordant TUB-band MTBDR_{sl} results on culture isolates (one TUB-positive and FQ and SLID sensitive, one TUB-positive and FQ and SLID resistant) where both CE results were TUB-band negative. There was also 1/20 (5%) discordant SLID result (CE showed resistance but the isolate showed susceptibility). Importantly, all three discordant results corresponded to a “very low” semiquantitation (C_T > 28.0). All TUB-band, susceptibility and resistance calls were concordant at C_T ≤ 24, indicating that the diagnostic accuracy of MTBDR_{sl} on CE vs. isolates is likely comparable at this threshold.

Spoligotyping on cartridge extract. *Dilution series.* Spoligotyping resulted in a readable strain type for dilutions ≥ 10³ CFU/ml, corresponding to the same threshold seen for MTBDR_{sl}.

Clinical specimens. Spoligotyping on specimen CE and crude DNA from matched culture isolates were highly concordant 10/10 (100%) at the threshold defined by the dilution series (Table 3). A variety of strain families were observed with Beijing as the predominant family type [6/10 (60%)] as well as 2/10 (20%) LAM and 2/10 (20%) T1 family type.

Targeted sequencing on extract from used Xpert MTB/RIF cartridges. *Dilution Series.* Targeted Sanger sequencing was done on dilution series CE. For the *rrs* PCR on CE, sequence shorter than the expected length was observed. PCR of *gyrA* from CE from dilutions 10³–10⁴ CFU/ml resulted in sequence expected length, however high background noise occurred and the sequence did not align to H37Rv [NC_000962]. *gyrA* on CE from dilutions 10⁵–10⁶ CFU/ml aligned to the reference genome, however, several SNPs known to be present in the resistance determining regions (identified by sequencing of the corresponding isolate) were not detected. Due to the relatively poor limit of detection and accuracy of Sanger sequencing on dilution series CE, we did not do sequencing on clinical specimen CEs.

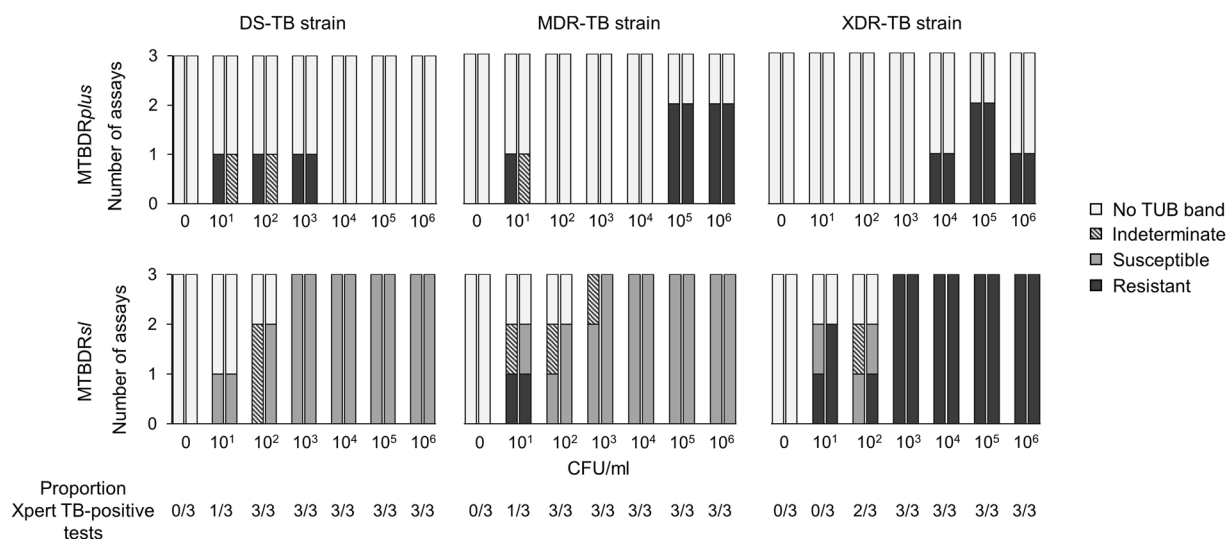


Figure 2. Results of MTBDR_{plus} and MTBDR_{sl} on Xpert CE from a dilution series of DS-, MDR- and XDR-TB strains. MTBDR_{plus} (irrespective of concentration and strain) had high TUB-band negativity and indeterminate rates. However, MTBDR_{sl} had high sensitivity and specificity and low indeterminate rates. For each dilution, left bars are for rifampicin (MTBDR_{plus}, top row) or fluoroquinolones (MTBDR_{sl}, bottom row) and right bars are for isoniazid (MTBDR_{plus}) or second-line injectables (MTBDR_{sl}). Data from LPA on DS-TB, MDR-TB and XDR-TB strains are shown. The experiment was done in triplicate. Abbreviations: CFU – colony forming; DS-TB – drug susceptible TB; MDR-TB – multidrug resistant TB; XDR-TB – extensively drug resistant TB; units; Xpert – Xpert MTB/RIF.

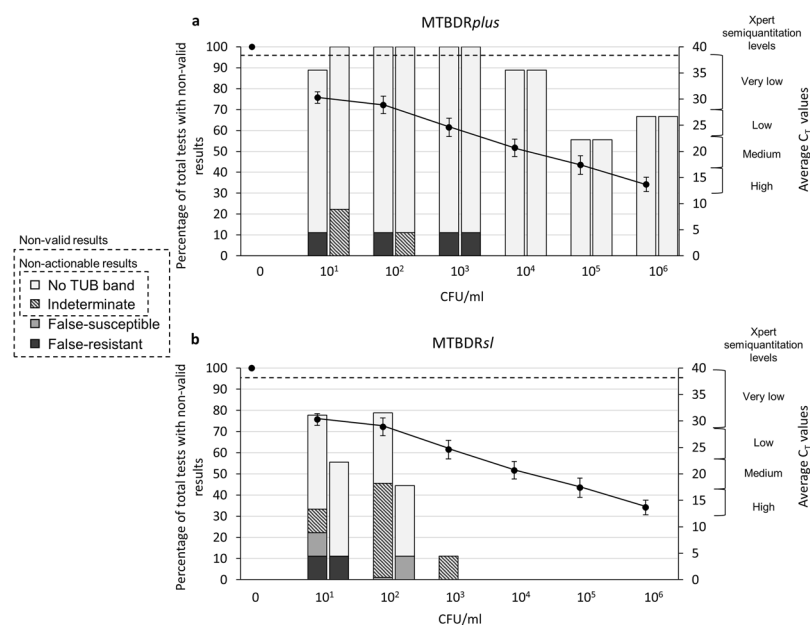


Figure 3. Xpert MTB/RIF quantitative information [average cycle threshold (C_T) values] (line graph, right y-axes) versus bacterial load (CFU/ml) in a triplicate dilution series for MTBDR_{plus} (a) and MTBDR_{sl} (b) done on CE. Left y-axes (bars) show the proportion of assays with non-valid results, disaggregated into non-actionable (TUB-band negative, indeterminate) and non-valid (false-susceptible, false-resistant). For each dilution, left bars are for rifampicin (MTBDR_{plus}, top) or fluoroquinolones (MTBDR_{sl}, bottom) and right bars are for isoniazid (MTBDR_{plus}) or second-line injectables (MTBDR_{sl}). Beyond 10³ CFU/ml, there were no false resistance or susceptibility calls for MTBDR_{sl}, which corresponds to $C_T \leq 24$. $C_T \geq 38$ (horizontal dashed line) correspond to a negative Xpert. Error bars show standard error (SE) of average C_T . Right y-axes show C_T corresponding to Xpert semiquantitation levels of very low ($C_T > 28$), low ($C_T = 22-28$), medium ($C_T = 16-22$) and high ($C_T < 16$). Pooled data from LPAs on DS-TB, MDR-TB and XDR-TB strains are shown. Abbreviations: CFU – colony forming units; DS-TB – drug susceptible TB; MDR-TB – multidrug resistant TB; XDR-TB – extensively drug resistant TB; CFU – colony forming units; Xpert – Xpert MTB/RIF.

| All Xpert TB-positive specimens | | | | | | Xpert TB-positive specimens with $C_T \leq 24$ | | | | | |
|-------------------------------------|-------------|-------------------------------------|----------------------------|-------------------------------------|------------|--|------------|--------------------------------------|----------------------------|--------------------------------------|------------|
| Xpert rifampicin-susceptible | | | Xpert rifampicin-resistant | | | Xpert rifampicin-susceptible | | | Xpert rifampicin-resistant | | |
| MTBDRplus (n = 56) | | MTBDRsl (n = 56) | | MTBDRsl (n = 29)* | | MTBDRplus (n = 49) | | MTBDRsl (n = 49) | | MTBDRsl (n = 20)* | |
| TUB-band positive (%) 47/56 (84) | | TUB-band positive (%) 47/54 (96) | | TUB-band positive (%) 27/29 (93) | | TUB-band positive (%) 45/49 (92) | | TUB-band positive (%) 49/49 (100) | | TUB-band positive (%) 20/20 (100) | |
| Rifampicin (%) | | Fluoroquinolones (%) | | Fluoroquinolones (%) | | Rifampicin (%) | | Fluoroquinolones (%) | | Fluoroquinolones (%) | |
| Susceptible | 0/47 (0) | Susceptible | 53/54 (98) | Susceptible | 14/27 (52) | Susceptible | 0/56 (0) | Susceptible | 49/49 (100) | Susceptible | 11/20 (55) |
| Resistant | 47/47 (100) | Resistant | 0/54 (0) | Resistant | 10/27 (37) | Resistant | 45/49 (92) | Resistant | 0/49 (0) | Resistant | 9/20 (45) |
| Indeterminate | 0/47 (0) | Indeterminate | 1/54 (2) | Indeterminate | 3/27 (11) | Indeterminate | 0/49 (0) | Indeterminate | 0/49 (0) | Indeterminate | 0/20 (0) |
| Isoniazid (%) | | Second-line injectables (%) | | Second-line injectables (%) | | Isoniazid (%) | | Second-line injectables (%) | | Second-line injectables (%) | |
| Susceptible | 11/47 (23) | Susceptible | 48/54 (88) | Susceptible | 15/27 (56) | Susceptible | 11/49 (23) | Susceptible | 46/49 (94) | Susceptible | 14/20 (70) |
| Resistant | 0/47 (0) | Resistant | 3/54 (6) | Resistant | 9/27 (33) | Resistant | 0/49 (0) | Resistant | 1/49 (2) | Resistant | 6/20 (30) |
| Indeterminate | 36/47 (77) | Indeterminate | 3/54 (6) | Indeterminate | 3/27 (11) | Indeterminate | 34/49 (69) | Indeterminate | 2/49 (4) | Indeterminate | 0/20 (0) |
| TUB-band negative (%) 9/56 (16) | | TUB-band negative (%) 2/56 (4) | | TUB-band negative (%) 2/29 (7) | | TUB-band negative (%) 4/49 (8) | | TUB-band negative (%) 0/49 (0) | | TUB-band negative (%) 0/20 (0) | |

Table 2. Results of MTBDRplus and MTBDRsl drug susceptibility testing using cartridge extract on clinical specimens. MTBDRplus had high indeterminate results and rifampicin-resistance false-positive rates. MTBDRsl had low indeterminate rates for both DS-TB and DR-TB specimens and performance improved when MTBDRsl was done only on specimens with $C_T \leq 24$. *For the 29 Xpert RIF-resistant specimens we were able to retrieve 20 paired culture isolates used for MTBDRsl. 18/20 matched for FQs and 17/20 for SLIDs, the 2/20 done on crude DNA had LPA results whereas the LPA on CE was TUB-band negative. 1/20 did not match for the SLID resistance. Both the TUB-band negative and discordant SLIDs result corresponded to “very low” semi-quantitation level. When defined threshold of $C_T \leq 24$ was applied all LPAs on CE matched LPA from culture isolates.

| Material used for spoligotyping | Xpert MTB/RIF semi-quantitation level | C_T | Spoligotyping pattern | Family |
|---------------------------------|---------------------------------------|-------|-----------------------|------------|
| Negative Control | | | N/A | N/A |
| H37Rv | | | | H37Rv |
| BCG | | | | BOVIS1_BCG |
| Isolate (a) | - | - | | T1 |
| CE (a) | MEDIUM | 20.85 | | T1 |
| Isolate (b) | - | - | | BEIJING |
| CE (b) | HIGH | 12.28 | | BEIJING |
| Isolate (c) | - | - | | BEIJING |
| CE (c) | HIGH | 16.20 | | BEIJING |
| Isolate (d) | - | - | | BEIJING |
| CE (d) | HIGH | 15.00 | | BEIJING |
| Isolate (e) | - | - | | BEIJING |
| CE (e) | LOW | 24.40 | | BEIJING |
| Isolate (f) | - | - | | BEIJING |
| CE (f) | HIGH | 14.74 | | BEIJING |
| Isolate (g) | - | - | | LAM9 |
| CE (g) | LOW | 20.86 | | LAM9 |
| Isolate (h) | - | - | | BEIJING |
| CE (h) | MEDIUM | 21.35 | | BEIJING |
| Isolate (i) | - | - | | LAM3 |
| CE (i) | HIGH | 16.28 | | LAM3 |
| Isolate (j) | - | - | | T1 |
| CE (j) | MEDIUM | 21.88 | | T1 |

Table 3. Spoligotyping results performed on CE done on sputum specimens and paired culture isolates at defined threshold ($C_T \leq 24$).

Discussion

Our key findings are: (1) MTBDR*sl* on CE enabled genotypic drug-susceptibility testing for FQs and SLIDs with high accuracy and low indeterminate rates when the Xpert semiquantitation category was at least “medium” or $C_T \leq 24$ (corresponding to $\geq 10^3$ CFU/ml), (2) spoligotyping was feasible and accurate on CE at the same threshold, (3) MTBDR*plus* was not feasible or accurate on CE and (4) neither was Sanger sequencing. These data have implications for the routine diagnosis of drug-resistant TB, researchers, and test developers.

Xpert is one of the most widely used tests for TB and drug-resistance^{9,34} and although it is a significant advancement, time-to-treatment – especially for MDR- and XDR-TB – is still very long^{35–38}. Our results show that accurate second-line drug testing using MTBDR*sl* is possible on CE from Xpert cartridges that would otherwise be discarded. This potentially allows for a rapid, single-specimen diagnosis of XDR-TB without additional specimen collection. Importantly, we defined a threshold at which this approach is feasible, meaning that MTBDR*sl* assays do not need to be wasted on CE unlikely to give a valid result. Using this threshold, we showed that on clinical specimen CEs, susceptibility and resistance calls were concordant with those from the isolate^{19,39}. Furthermore, we showed that it is possible to do spoligotyping on CE at this threshold, which will inform strain surveillance and research studies on relapse and reinfection where specimens are limited. Collectively, these findings may reduce the need for culture.

Although our data suggest that the MTBDR*sl* will work on CE from cartridges with an Xpert semiquantitation category of at least “low”, we suggest that, in laboratories where C_T cannot be readily calculated, a category of at least “medium” is used to guide use of this strategy unless the laboratory is comfortable with some semiquantitation low specimens not having a valid MTBDR*sl* result. Alternatively, if smear microscopy is available, smear-positivity may be used to guide use of CE, however, some smear-negative specimens in whom this approach would work (10^3 – 10^4 CFU/ml) would be unnecessarily excluded.

When considering the CE approach, it is important to identify a safe and sterile environment to avoid contamination. Although Xpert sample reagent as well as the sonication lysis step within the cartridge helps ensure *M. tuberculosis* is no longer culturable (and therefore poses minimal infectious risk⁴⁰), steps to minimise the risk of *rpoB* amplicon cross-contamination should be implemented. These can include working in a dedicated cabinet or room and sterilising the work area with UV and disinfectant after CE is collected. Importantly, however, cross-contamination of other Xpert cartridges with *rpoB* amplicons appears unlikely. Although Xpert's automated pre-amplification wash step does not remove large pieces of debris-associated genomic DNA, it does efficiently remove high concentrations of contaminating *rpoB* amplicons from assays like MTBDR*plus*^{41,42}. NAATs without such a wash step may be more vulnerable to CE *rpoB* amplicon cross-contamination.

Our study differed from a previous study which showed that sequencing, MTBDR*plus*, spoligotyping and MIRU-VNTR typing are feasible on the sputum mixed with Xpert sample reagent⁴³. However, this sample reagent method has a number of disadvantages: 1) often no volume remains, 2) prolonged exposure to sample reagent degrades DNA and potentially introduces mutations^{9,40}, and 3) it still requires DNA extraction prior to PCR. Furthermore, DNA extraction adds cost and is not always feasible in laboratories in high burden countries; whereas the CE method yields directly usable material and does not need additional extraction or purification steps. An advantage, however, of using the sputum mixed with Xpert reagent buffer, is that it likely avoids high MTBDR*plus* error rates (TUB-band negative, indeterminate, false-positive) seen with CE. This could be due to the large amount of *rpoB* amplicons in Xpert TB-positive CE, which share binding sites with MTBDR*plus* probes and confound the assay resulting in non-valid results. Furthermore, the *rpoB* PCR that occurs as part of MTBDR*plus* may sequester reagents away from the multiplex *inhA* and *katG* amplification reactions. Testing for mutations conferring INH resistance using CE might hence be possible with the Genoscholar INH II line probe assay (which does not contain *rpoB* probes)⁴⁴. Sequencing from CE thus primarily appears to be driven by *rpoB* amplicon interference (although a PCR clean-up was done prior to sequencing, this would have co-purified *rpoB* amplicons). Further investigation with primers optimised for minimal-input DNA may be warranted, however, it appears that, for sequencing, the best approach to avoid contaminating amplicons might be to PCR from the specimen-Xpert sample reagent mixture⁴⁵. Given the rates of non-valid CE results below the defined threshold, we suggest that specimen-Xpert sample reagent mix be kept in the event that C_T falls >24 .

The results presented here should be interpreted in context of their limitations. For the clinical specimens tested from the NHLS, matched culture isolates were not available for Xpert RIF-susceptible specimens, as per the national algorithm. However, the dilution series experiments showed very high concordance between MTBDR*sl* on CE vs. the isolates. The utility of CE depends on the downstream test used and MTBDR*sl* susceptible or non-valid results should be interpreted from CE the same as when they are done on patient specimens (i.e., further investigation, including culture, is recommended)⁴⁵. Realistically, cartridges may need to be transported from remote locations and so the effect of storing cartridges for prolonged duration (>5 days) and at ambient temperature requires further systematic testing. Using bacilli in buffer can have limitations, which is why we also used patient clinical specimens, which are a better material to test than bacilli added to sputum (the former has bacilli within a sputum matrix, whereas in the latter bacilli are typically freely floating in bubbles).

In conclusion, CE contains template DNA for second-line DST using MTBDR*sl*, resulting in accurate results highly concordant with those from isolates, provided bacillary load in the specimen corresponds to at least a “medium” Xpert semiquantitation category of $C_T \leq 24$. This potentially facilitates XDR-TB detection within days from a single specimen. Spoligotyping is also feasible on CE and works consistently at this threshold. Our method provides an opportunity to potentially reduce the burden associated with additional specimen collection, such as patient treatment delay, pre-treatment loss-to-follow-up, and increased patient and provider costs. Furthermore, it shows that material that would otherwise be discarded still holds diagnostic utility.

References

- World Health Organization. Global tuberculosis report (2016).
- Pooran, A., Pieterse, E., Davids, M., Theron, G. & Dheda, K. What is the cost of diagnosis and management of drug resistant tuberculosis in South Africa? *PloS one* **8**, e54587 (2013).
- Cox, H., Ramma, L., Wilkinson, L., Azevedo, V. & Sinanovic, E. Cost per patient of treatment for rifampicin-resistant tuberculosis in a community-based programme in Khayelitsha, South Africa. *Tropical medicine & international health: TM & IH* **20**, 1337–1345, <https://doi.org/10.1111/tmi.12544> (2015).
- Dheda, K. *et al.* The epidemiology, pathogenesis, transmission, diagnosis, and management of multidrug-resistant, extensively drug-resistant, and incurable tuberculosis. *The Lancet Respiratory Medicine* **5**, 291–360 (2017).
- Cox, H. *et al.* Community-based treatment of drug-resistant tuberculosis in Khayelitsha, South Africa. *The International Journal of Tuberculosis and Lung Disease* **18**, 441–448 (2014).
- Moyo, S. *et al.* Loss from treatment for drug resistant tuberculosis: risk factors and patient outcomes in a community-based program in Khayelitsha, South Africa. *PloS one* **10**, e0118919 (2015).
- US Food Drug Administration. New data shows test can help physicians remove patients with suspected TB from isolation earlier (2015).
- Boehme, C. C. *et al.* Rapid molecular detection of tuberculosis and rifampin resistance. *New England Journal of Medicine* **363**, 1005–1015 (2010).
- World Health Organization. Xpert MTB/RIF implementation manual: technical and operational 'how-to': practical considerations [Internet]. Geneva: World Health Organization (2014).
- World Health Organization. *Xpert Rollout Status* (2016). accessed at: <http://www.who.int/tb/areas-of-work/laboratory/mtb-rif-rollout/en/>, 19 June 2017.
- Health Department of the Republic of South Africa. National Tuberculosis Management Guidelines (2014).
- Sachdeva, K. S. Management of Tuberculosis: Indian Guidelines. *API Medicine Update, Section 15*, 479e483 (2013).
- Jacobson, K. R. *et al.* Implementation of GenoType MTBDR plus Reduces Time to Multidrug-Resistant Tuberculosis Therapy Initiation in South Africa. *Clinical infectious diseases* **56**, 503–508 (2012).
- MacPherson, P., Houben, R. M., Glynn, J. R., Corbett, E. L. & Kranzer, K. Pre-treatment loss to follow-up in tuberculosis patients in low-and lower-middle-income countries and high-burden countries: a systematic review and meta-analysis. *Bulletin of the World Health Organization* **92**, 126–138 (2014).
- Peter, J. G. *et al.* Comparison of two methods for acquisition of sputum samples for diagnosis of suspected tuberculosis in smear-negative or sputum-scarce people: a randomised controlled trial. *The Lancet Respiratory Medicine* **1**, 471–478 (2013).
- Sharath, B. & Shastri, S. India's new TB diagnostic algorithm-far from reality? *Public health action* **6**, 206–206 (2016).
- Peter, J. G., Theron, G., Singh, N., Singh, A. & Dheda, K. Sputum induction to aid diagnosis of smear-negative or sputum-scarce tuberculosis in adults in HIV-endemic settings. *European Respiratory Journal* **43**, 185–194 (2014).
- Tomasichio, M. *et al.* The diagnostic accuracy of the MTBDRplus and MTBDRsl assays for drug-resistant TB detection when performed on sputum and culture isolates. *Scientific reports* **6** (2016).
- Theron, G. *et al.* The diagnostic accuracy of the GenoType () MTBDRsl assay for the detection of resistance to second-line anti-tuberculosis drugs. *Cochrane Database of Systematic Reviews* **10** (2014).
- Nathavitharana, R. R. *et al.* Accuracy of line probe assays for the diagnosis of pulmonary and multidrug-resistant tuberculosis: a systematic review and meta-analysis. *European Respiratory Journal* **49**, 1601075 (2017).
- Global Laboratory Initiative. Practical Guide to TB Laboratory Strengthening (2017).
- Metcalfe, J. Z. *et al.* Mycobacterium tuberculosis subculture results in loss of potentially clinically relevant heteroresistance. *Antimicrobial agents and chemotherapy*, AAC. 00888–00817 (2017).
- Dheda, K. *et al.* TB drug resistance in high-incidence countries. *Tuberculosis* **58**, 95–110 (2012).
- Dlamini-Mvelase, N. R., Werner, L., Philo, R., Cele, L. P. & Mlisana, K. P. Effects of introducing Xpert MTB/RIF test on multi-drug resistant tuberculosis diagnosis in KwaZulu-Natal South Africa. *BMC infectious diseases* **14**, 442 (2014).
- Cox, H. *et al.* Delays and loss to follow-up before treatment of drug-resistant tuberculosis following implementation of Xpert MTB/RIF in South Africa: A retrospective cohort study. *PLoS medicine* **14**, e1002238 (2017).
- Engel, N. *et al.* Compounding diagnostic delays: a qualitative study of point-of-care testing in South Africa. *Tropical medicine & international health: TM & IH* **20**, 493–500, <https://doi.org/10.1111/tmi.12450> (2015).
- Naidoo, P. *et al.* Estimation of losses in the tuberculosis care cascade in South Africa and methodological challenges. *Under revision*.
- Truant, J., Brett, W. & Thomas, W. Jr Fluorescence microscopy of tubercle bacilli stained with auramine and rhodamine. *Henry Ford Hospital Medical Bulletin* **10**, 287–296 (1962).
- Hain Lifescience. GenoType MTBDRplus VER 2.0 Instructions for Use (2014).
- Hain Lifescience. GenoType MTBDRsl VER 2.0 Instructions for Use (2015).
- Streicher, E. *et al.* Spoligotype signatures in the Mycobacterium tuberculosis complex. *Journal of clinical microbiology* **45**, 237–240 (2007).
- Kamerbeek, J. *et al.* Simultaneous detection and strain differentiation of Mycobacterium tuberculosis for diagnosis and epidemiology. *Journal of clinical microbiology* **35**, 907–914 (1997).
- Demay, C. *et al.* SITVITWEB—a publicly available international multimer database for studying Mycobacterium tuberculosis genetic diversity and molecular epidemiology. *Infection, Genetics and Evolution* **12**, 755–766 (2012).
- Lawn, S. D. & Nicol, M. P. Xpert® MTB/RIF assay: development, evaluation and implementation of a new rapid molecular diagnostic for tuberculosis and rifampicin resistance. *Future microbiology* **6**, 1067–1082 (2011).
- Metcalfe, J. Z. *et al.* Xpert® MTB/RIF detection of rifampin resistance and time to treatment initiation in Harare, Zimbabwe. *The International Journal of Tuberculosis and Lung Disease* **20**, 882–889 (2016).
- Iruedo, J., O'Mahony, D., Mabunda, S., Wright, G. & Cawc, B. The effect of the Xpert MTB/RIF test on the time to MDR-TB treatment initiation in a rural setting: a cohort study in South Africa's Eastern Cape Province. *BMC infectious diseases* **17**, 91 (2017).
- Cox, H. S. *et al.* Impact of Xpert MTB/RIF for TB diagnosis in a primary care clinic with high TB and HIV prevalence in South Africa: a pragmatic randomised trial. *PLoS medicine* **11**, e1001760 (2014).
- van Kampen, S. C. *et al.* Effects of introducing Xpert MTB/RIF on diagnosis and treatment of drug-resistant tuberculosis patients in Indonesia: a pre-post intervention study. *PloS one* **10**, e0123536 (2015).
- Theron, G. *et al.* GenoType® MTBDRsl assay for resistance to second-line anti-tuberculosis drugs. *The Cochrane Library* (2016).
- Banada, P. P. *et al.* Containment of bioaerosol infection risk by the Xpert MTB/RIF assay and its applicability to point-of-care settings. *Journal of clinical microbiology* **48**, 3551–3557 (2010).
- Theron, G. *et al.* Xpert MTB/RIF results in patients with previous tuberculosis: can we distinguish true from false positive results? *Clinical Infectious Diseases* **62**, 995–1001 (2016).
- Blakemore, R. *et al.* Evaluation of the analytical performance of the Xpert MTB/RIF assay. *Journal of clinical microbiology* **48**, 2495–2501 (2010).
- Alame-Emane, A. K. *et al.* The use of GeneXpert remnants for drug resistance profiling and molecular epidemiology of tuberculosis in Libreville, Gabon. *Journal of clinical microbiology*, JCM. 02257–02216 (2017).
- Mitarai, S. *et al.* Comprehensive multicenter evaluation of a new line probe assay kit for identification of Mycobacterium species and detection of drug-resistant Mycobacterium tuberculosis. *Journal of clinical microbiology* **50**, 884–890 (2012).
- World Health Organization. The use of molecular line probe assays for the detection of resistance to second-line anti-tuberculosis drugs: policy guidance (2016).

Acknowledgements

The authors would like to thank Dr Michael Whitfield who provided the strains for the *in vitro* study. The financial assistance of the National Research Foundation (NRF) towards this research is hereby acknowledged. Opinions expressed and conclusions arrived at, are those of the author(s) and are not necessarily to be attributed to the NRF. Research reported in this publication was supported by the South African Medical Research Council. The content is solely the responsibility of the authors and does not necessarily represent the official views of the South African Medical Research Council. Research reported in this publication was also supported by the Faculty of Medicine and Health Sciences, Stellenbosch University and the European & Developing Countries Clinical Trials Partnership 2.

Author Contributions

G.T., R.W. and M.D.V. conceived the experiments. R.V., B.D., S.P., N.K. and A.R. conducted the experiments. J.S. and T.D. provided specimens and data from the NHLS. R.V. and B.D. analysed the data. All authors reviewed the manuscript.

Additional Information

Supplementary information accompanies this paper at <https://doi.org/10.1038/s41598-017-14385-x>.

Competing Interests: The authors declare that they have no competing interests.

Publisher's note: Springer Nature remains neutral with regard to jurisdictional claims in published maps and institutional affiliations.



Open Access This article is licensed under a Creative Commons Attribution 4.0 International License, which permits use, sharing, adaptation, distribution and reproduction in any medium or format, as long as you give appropriate credit to the original author(s) and the source, provide a link to the Creative Commons license, and indicate if changes were made. The images or other third party material in this article are included in the article's Creative Commons license, unless indicated otherwise in a credit line to the material. If material is not included in the article's Creative Commons license and your intended use is not permitted by statutory regulation or exceeds the permitted use, you will need to obtain permission directly from the copyright holder. To view a copy of this license, visit <http://creativecommons.org/licenses/by/4.0/>.

© The Author(s) 2017